

Elite

Transport Wheelchair

TRANSPORT
BRODA[®]



Shown with armrests, legrests, rear mag wheels, and solid seat pan.

Patient, Facility, & Caregiver Benefits

- Facilitates positioning for comfort and respiratory function.
- Tilt and recline provides ease of use for frequent re-positioning.
- Proper seating alignment decreases postural deviations.

FOR MORE INFORMATION, VISIT BRODASEATING.COM




Rear mount mag wheels facilitate hand propulsion allowing patients to engage in their surroundings



Adjustable tilt and recline helps decrease postural deviations such as slumping, lateral lean and head drop



Posterior tilt effectively opens the diaphragm improving digestion, oxygenation, and organ function

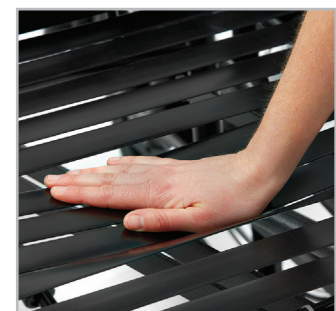


Comfort Tension Seating[®] improves postural support and sitting tolerance

VERSION 1.0 05/2020

UPHOLSTERY OPTIONS



FRAME OPTIONS

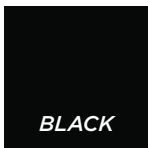

Industry-leading Positioning, Versatility, and Function.

SPECIFICATIONS

Degrees of Tilt	up to 40°
Degrees of Recline	up to 90°
Weight Capacity	350 lbs. (158 kg)
Seat Widths	16", 18", 20", 22", 24", 26", 28" (40cm, 45 cm, 50cm, 55cm, 60cm, 65cm)
Seat Depths	15.5", 17", 18.5", 20" (37.5cm, 43cm, 46cm, 50cm)
Floor-to-Seat Heights	19", 20.5", 22", 23.5" (48cm, 51cm, 55cm, 59cm)
Wheels	5" (12cm) Total Lock Front Casters 20" or 22" (50cm or 55cm) Rear Mags
Lifetime Warranty Frame	10 years on 16 Gauge tubular steel frame
Pressure-Reducing Seat & Backrest	Integrated with Comfort Tension Seating®
Standard Padding*	Lower side pads, upper side pads, seat and back pad with cover

Additional accessories available: foot supports, specialty padding packages, trays, and positioning straps